Professional Sonic Toothbrush

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Professional Sonic Toothbrush

Instructions for use and operation

Please read this information carefully before using the product and keep it for further reference.

© Copyright: ApaCare is a trademark of Cumdente GmbH. All rights reserved.

Status: February 2022

Important safety information

- Please do not charge near water-bearing devices.
- The device contains a built-in rechargeable battery.
- Avoid heat, open flames and magnetic fields.
- Check device before use.
- Please do not use defective devices.
- Do not open the device.
- Have repairs carried out only by the manufacturer.
- Exercise caution after dental-surgical therapy. Wait for 3 months and then consult your dentist before use.
- First-time users may experience discomfort due to the rapid brush oscillation. This discomfort should disappear after no more than 2 weeks.
- Children under 8 years of age, patients in need of care or disabled patients should not use the device without supervision or use by parents/carers.

Please keep out of reach of children. Toothbrushes are patientspecific devices. Please do not use them together.

Specifications

ApaCare[®] Sonic Toothbrush


USB Charger

Voltage: 5V
Power: 1W
Output: 5V, 55mA

Device Voltage: 3,7V
Power: 1,2W
Standard: GB4706.1_2005 GB4706.59-2008
Country of Origin: People's Republic of China

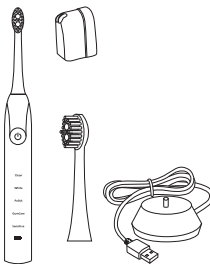Scope of delivery

- 1x Handpiece
- 2x Brush head with protective cover
- 1x Charging station with USB cable and power adapter 5V, 1W

Functional description (Mode)

-  **Clean:** Usual daily cleaning
-  **White:** Whitening of teeth, especially for smokers or when drinking tea or coffee
-  **Polish:** Polishing of the tooth surface
-  **Gumcare:** Gentle gum care
-  **Sensitive:** Gentle mode, especially for sensitive patients or children

Operation

Switch on

Press the switch-on button once. The toothbrush starts the set mode.


Switch off

Switching off by pressing the function key again after at least 5 seconds of operation. The unit switches off.


or

Press the function key (Mode) once more after the 5th press. The unit switches off.


Directions for use

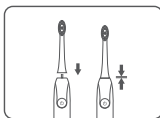
1. **Attach the brush head** to the toothbrush by slipping it onto the shaft (see illustration). A gap of 0.5 mm - 1 mm remains between the brush head and the toothbrush housing, which is required for unhindered oscillation of the brush head.

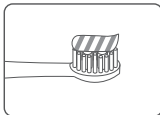
2. **Apply a small amount of toothpaste** (e.g. ApaCare[®] toothpaste) onto the brush head.

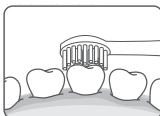
3. **Place the brush with the brush head against the row of teeth** and switch on the toothbrush.

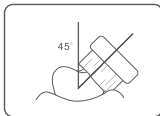4. **Place the brush head at a 45° angle** to the teeth and gums and press gently. Make small circular or oscillating movements. Do not scrub. Avoid high contact pressure. Carefully clean all the teeth of each jaw half in groups of approx. 3 teeth inside, outside and on the occlusal surface (After at least 30s per jaw half, „Interval-Timer“ will set in).

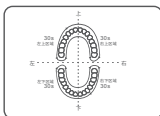5. **After use, rinse the toothbrush** thoroughly with water and store to dry (without brush cover). Clean toothbrush shaft regularly after brush is removed.

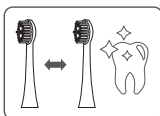6. **Brush your teeth in this way at least 2 times a day.**


7. **Replace the brush head at least every 3 months.** If the bristles become deformed quickly, be careful not to use excessive contact pressure hence avoid premature replacement of the brush head.

8. **Each brush head is intended for individual personal use.** It must not be used by more than one person. If necessary, label brush heads to personalize them.

Charging


1. Fully charge your new toothbrush by plugging it onto the USB charger (12 hours) (full charge is indicated by the green indicator light on the charger becoming permanently lit). Do this also if you have not used the toothbrush in more than 3 months or if the vibration seems insufficient.
2. Place the charging station in a safe, level, dry and cool place and make sure that it cannot fall into water or onto heat-bearing surfaces.
3. Plug in the USB charging cable.
4. Place the toothbrush vertically on the charging unit. Switch off the toothbrush for charging. The indicator light on the charger flashes during charging.
5. Do not switch on the toothbrush during charging.
6. When charging is complete, it is automatically terminated.
7. A flashing indicator light on the toothbrush shows that the toothbrush must be recharged.

Cleaning

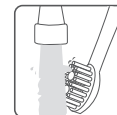1. After each use, rinse the brush head thoroughly under running water and allow to air dry.


2. Do not pull on the brushes.

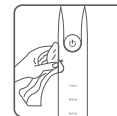3. Remove the brush head from the toothbrush at least once a week and clean the connection points carefully.

4. Clean the toothbrush with running water and dry, for example with a soft towel.

5. Do not immerse the toothbrush in water or cleaning agents.

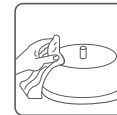6. Do not use hot water (>40°C), solvents or alkaline cleaning agents for cleaning.

7. Wipe the charging station with a soft cotton cloth and clean the metal contacts with a dry cotton cloth every 6 months.

Cumdente GmbH
Paul-Ehrlich-Strasse 11
72076 Tübingen, Germany

Manufacturer: Ningbo Xoralcare Electronic Technology Co., LTD
No. 1, Yangming Street Yuyao, Zhejiang, China

Troubleshooting

Error	Cause	Troubleshooting	Comments
Reduced vibration intensity	<ol style="list-style-type: none"> 1. Toothbrush pressed down too hard. 2. Toothbrush attachment is too tight (touches the toothbrush housing). 	<ol style="list-style-type: none"> 1. Reduce contact pressure. 2. Remove brush head from shaft and carefully reattach. Make sure there is a gap of 0.5 to 1 mm. 	<p>For questions please contact: ApaCare® Customer Service mail@apacare.de Or call Phone: +49 70 71 97 55 721</p>
Premature switch-off after a short time	<ol style="list-style-type: none"> 1. Charging status insufficient or unit not used for a long period of time. 2. Device not placed correctly (vertically) on charging station. 3. Defective USB charging adapter. 	<ol style="list-style-type: none"> 1. Charge the unit for 12 hours. 2. Place the unit vertically on the charging station. Observe flashing charging indicator. 3. Battery life has expired. Dispose of device in accordance with regulations. 	
Cleaning process is interrupted	<ol style="list-style-type: none"> 1. Battery dead. 2. Malfunction. 	<ol style="list-style-type: none"> 1. Charge device. 2. Contact ApaCare® customer service. 	

Thank you for choosing our product.
 We will provide any support you need with our After Sales Service.

Warranty

The warranty period for the ApaCare® Sonic toothbrush is 1 year from the date of purchase. Brush heads are generally excluded from the warranty.

In the event of justified complaints during the warranty period, the toothbrush will be replaced.

The warranty is excluded in the following cases:

1. Incorrect use of the toothbrush.
2. Damage after attempting to open the device by yourself or attempting to repair the device by yourself.
3. Self-inflicted damage.
4. Damages related to force majeure (fire, natural disasters, ...).
5. In case of signs of wear or scratches on the surface.
6. Damage due to improper use, e.g. industrial use.
7. Repairs or service after the warranty period will be charged.
8. For all questions please contact the ApaCare® customer service: mail@apacare.com or Tel.: +49 70 71 97 55 721

To the extent that exclusion of liability is permitted by law, the manufacturer and Cumdente GmbH shall have no liability whatsoever for any loss or damage caused by this product, whether direct, indirect, special, incidental or consequential, regardless of the legal basis, including warranty, contract, negligence or intent.

Note!

Replace brush heads every 3 months at the latest.

ApaCare® Sonic 3D Sensitive Brushes
 2 pieces

Art.-No.: 1002022

Available at: www.apacare.com

ApaCare® Home

Dental care products with liquid enamel

By using ApaCare® dental care products, you are following the recommendation of many dentists. All ApaCare® dental care products contain liquid tooth enamel of similar composition to natural teeth (medical hydroxyapatite). For naturally smooth, healthy and white teeth.

CLEAN

REPAIR

PROTECT

WHITEN

ApaCare® Toothpaste

Remineralizing toothpaste with medical hydroxyapatite and fluoride

- ✓ 2-3 times daily

ApaCare® Repair

Agent Repair

- ✓ Daily after brushing teeth

ApaCare® Polish

Tooth polishing paste

- ✓ 1-2 times / week

ApaCare® Liquid

Antibacterial tooth/mouth rinse with medical hydroxyapatite

- ✓ 2 times a day, 30 s

Warranty card

(Please send to ApaCare® customer service)

Model no: _____

Date: _____

Name: _____

Tel. no: _____

Address: _____

Error description: _____

Cumdente GmbH
 ApaCare® Customer Service
 Paul-Ehrlich-Strasse 11
 72076 Tübingen, Germany

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© Copyright: ApaCare ist eine Marke der Cumdente GmbH. All rights reserved.